

ASSISTIVE TECHNOLOGY  
PRODUCT CATALOG  
2025-2026



## QUHA ZONO X - GYROSCOPIC BLUETOOTH® HEAD MOUSE



Quha Zono X is a tiny handsfree Bluetooth® head mouse for people who can't use desktop mice with their hands. With simple head movements the user can move the cursor accurately on any modern computer or mobile device.

Quha Zono X users can perform all regular mouse functions without hands. Quha Zono X is compatible with all modern mouse operated devices that support Bluetooth®: computers, tablets and smartphones.

Gesture recognition reads and interprets head movements to activate common features such as cursor pausing, window scrolling, cursor centering and multi-device switching.

The vertical and horizontal cursor speed can be adjusted separately while unwanted tremors can be reduced with the built-in tremor filter. The settings application is available for Windows and Android devices.

Quha Zono X has use time up to 18 hours with an easy-to-use fast charging dock. It also includes an audio indicator for confirming user input.

Sales package includes an easy-to-use magnetic ear clip and eyeglass clip for wearing the mouse.

Quha Zono X is compatible with all Bluetooth® enabled computers and mobile devices. Pair with several devices and switch between them with a gesture. Compatible operating systems include Windows, macOS, ChromeOS, Android, iPadOS, iOS and Linux.



## Quha Zono X mouse technical specifications

- Dimensions: Length 44 mm / 1.7 inches, width 25 mm / 1.0 inches, thickness 13 mm / 0.5 inches, weight 12 g / 0.4 ounces
- Bluetooth range: Up to 10 meters / 30 ft
- Battery: LiPo (Lithium Polymer) 90 mAh. Active usage time is up to 18 hours depending on power saving/ performance mode and use environment.
- Mouse refresh rate: 208 Hz

## Sales package includes:

- Quha Zono X mouse
- Charging dock
- License for Quha Dwell 2 software (page 10)
- USB-C charging cable
- Ear clip
- 3 silicone ear tips
- Eyeglass clip
- User manual and warranty card

Quha Zono X is EU MDR and UKCA compliant.

See additional ways to wear Quha Zono X on page 4.



P/N	PRODUCT
23	QUHA ZONO X GYROSCOPIC BLUETOOTH® HEAD MOUSE

# HEADBAND FOR QUHA ZONO X

Headband for Quha Zono X is a convenient wearing option for those who prefer the headband attachment and/or want to connect a switch to their Zono X mouse for clicking.

The headband features a magnetic holder for the Zono X mouse for quick and easy attachment. The holder also includes a 3.5 mm connector for connecting a switch for performing the functions of the left mouse button.

Quha Zono X can be used either on the left or right side.

P/N	PRODUCT
173	HEADBAND FOR QUHA ZONO X



# QUHA SENTO FOR ZONO X WITH ADJUSTABLE HEADBAND

Quha Sento for Zono X is a contactless puff switch. It is used to perform a short mouse click or dragging together with the Quha Zono X mouse. Mouse clicks are made easily and hygienically with a small puff to the microphone type of mouthpiece.

Using Sento is easy to learn and adapt to your personal puff. It is an essential addition to your Quha Zono X.

P/N	PRODUCT
195	QUHA SENTO FOR ZONO X WITH ADJUSTABLE HEADBAND





# QUHA ZONO 2 - ENHANCED GYROSCOPIC WIRELESS MOUSE

Quha Zono 2 is a wireless handsfree mouse operated with head movements. It detects motion and transforms it into intuitive onscreen cursor movement. Quha Zono 2 is an improved version of the Quha Zono mouse including all its features and a number of new ones.

Quha Zono 2 is for anyone who has difficulties in operating a computer mouse with hands. Typically Zono 2 users are people with various medical conditions relating to hand/arm functionality.

QUHA ZONO 2 INCLUDES INTEGRATED CURSOR CENTERING, PAUSING AND SCROLLING FEATURES.



Sales package includes:

- Quha Zono 2 mouse
- USB Receiver
- Head Mount Kit (see page 7)
- License for Quha Dwell 2 (see page 10)
- USB charging cable
- Stereo button adapter for Zono
- Adapter to connect two switches
- User manual and warranty card

Quha Zono 2 is compatible with all current Zono accessories. Choose the best way to wear Quha Zono 2 from the options presented on page 7.



//  
Love to have a single receiver with  
which you can choose if the  
switch is click, center or scroll!

P/N	PRODUCT
21	QUHA ZONO 2 GYROSCOPIC MOUSE

# QUHA ZONO - GYROSCOPIC WIRELESS MOUSE



Quha Zono makes it possible to use a computer in a flexible way, when you cannot use an ordinary mouse. It is an extremely light device, and since it is completely wireless there are endless ways to wear it.

Quha Zono is a gyrosopic mouse designed for anyone who has difficulties in operating a computer mouse with hands. It provides accurate and intuitive computer access with just small head movements.

WITH QUHA ZONO HEAD MOUSE YOU CAN DO EVERYTHING YOU CAN DO WITH THE MOST ACCURATE DESKTOP MICE.



Quha Zono is compatible with all Windows, macOS, Chromebook and Linux computers, iPads, Android tablets and most smartphones.

Sales package includes:

- Quha Zono mouse
- USB Receiver
- USB charging cable
- Stereo button adapter for Zono
- Adapter to connect two switches
- User manual and warranty card

Choose the best way to wear Quha Zono from the options presented on page 7.

//

I have been waiting for this mouse for 10 years.

P/N	PRODUCT
20	QUHA ZONO GYROSCOPIC MOUSE

# WAYS TO WEAR FOR QUHA ZONO AND ZONO 2

## HEADBAND

Headband is an easy way to wear Quha Zono and Zono 2. It is light and fits comfortably for hours. Headband is adjustable to make it suitable for most users.

Headband is an ideal alternative when using Zono while resting the head on a pillow or a headrest.

The Quha Zono mouse can be charged while attached to Headband.



## EYEWEAR KIT

Eyewear kit is a spectacle type way to wear Quha Zono and Zono 2.



## UNIVERSAL WEAR KIT

Universal Wear Kit enables Quha Zono and Zono 2 usage by hand or by foot. It includes a Velcro band and a soft holder for Quha Zono.



## HEAD MOUNT KIT

Head Mount Kit fits comfortably on most users' heads and can be worn for hours.



## BASEBALL CAP

Quha Baseball cap includes a classic cap and a white soft holder for Quha Zono and Zono 2.



## LIGHT EYEGLASS CLIP

Light Eyeglass clip is an ultra light clip to attach Quha Zono and Zono 2 safely to your own eyeglasses.



P/N	PRODUCT
170	HEADBAND, ADJUSTABLE
101	HEAD MOUNT KIT
103	LIGHT EYEGLASS CLIP
130	EYEWEAR KIT
111	UNIVERSAL WEAR KIT
132	BASEBALL CAP

## QUHA ZONO USER STORIES

### "IT IS VERY IMPORTANT THAT I CAN USE A COMPUTER AND A CAMERA BY MYSELF"

My name is Jari Mönkkönen and I am tetraplegic. I worked as a firefighter for 12 years. I also have a long career in artistic gymnastics. I have 29 gold medals at the Finnish championships and also a silver medal at the European championship on high-bar in '98.

My tetraplegia is caused by a sports accident. I made a double backflip and it landed on my neck. Nowadays I work as a head coach in a company called Voltimedia. We train young people and sports clubs.

Computers are very important in my daily life. After this accident, in a rehabilitation center, I tried several different devices and chose the Quha Zono head mouse. Quha Zono is accurate, fast, light-weight and easy to use.



### "I WOULD DEFINITELY RECOMMEND THIS DEVICE TO OTHERS"



My name is Katherine Quinn. I am 18 years old and I have cerebral palsy. My muscles don't work right, but my brain works just fine. I live in Carmichael, California. My interests are scrapbooking, horseback riding and watching videos. I'm going to be a senior at Rio Americano High school, where my favorite class is music. One day I would like to be Judge Katherine Quinn.


I have previously tried a headmouse, but it was difficult to use. I also tried the Eye Max, but it was too hard to set up every time I changed classes. Then I heard about the Quha products from Maggie Mahoney at Quha USA.

With the Quha Zono mouse I can easily upload a video on YouTube and read ebooks all by myself. Another great thing about the Quha Zono mouse is that you can move around with it – this is very important for me as sometimes my body moves by itself and I can't control it. These are the reasons why Zono is so great. I would definitely recommend this device to others.

SEE ALL USER STORIES AT [www.quha.com/user-stories](http://www.quha.com/user-stories)



# QUHA SENTO - CONTACTLESS PUFF SWITCH FOR QUHA ZONO AND ZONO 2

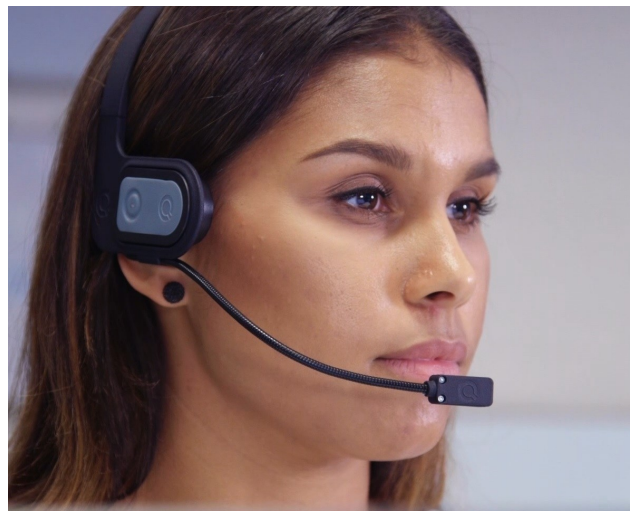

//

The design and the feeling of quality are great!

Quha Sento is the world’s first contactless puff switch. It is used to perform a short mouse click or dragging together with the Quha Zono and Zono 2 gyroscopic mice.

Mouse clicks are made easily and hygienically with a small puff to the microphone type of mouthpiece.

WITH QUHA SENTO YOU CAN PERFORM MOUSE CLICKS EASILY, HYGIENICALLY AND WITHOUT HANDS.

- DESIGNED CONSUMER ELECTRONICS LOOK & FEEL

  - No internal mouth pieces
- THE ONLY 100% HYGIENIC PUFF SWITCH ON THE MARKET

  - Enables one device to be shared between multiple users.
- COMPLETE WIRELESS COMPUTER ACCESS

  - Sento together with Quha Zono is the complete solution for computer access.
- NO CONTROL DELAY, UNIQUE SENSOR TECHNOLOGY

  - Instant puff response

P/N	PRODUCT
190	QUHA SENTO WITH ADJUSTABLE HEADBAND
191	QUHA SENTO WITH HEAD MOUNT KIT
192	QUHA SENTO WITH EYEWEAR KIT


# QUHA DWELL 2 - DWELL CONTROL FOR COMPLETE COMPUTER ACCESS

Quha Dwell 2 is a complete computer access software application for assistive mice.

The unique circular Instant Dwell feature makes computer usage faster, more comfortable and ergonomic compared to any current dwell control software.

With Quha Dwell 2 you can easily control all the functions of your computer, switch between applications, copy-paste and scroll. New functions can be added to Quha Dwell 2 to meet your needs.

P/N	PRODUCT
270	QUHA DWELL 2

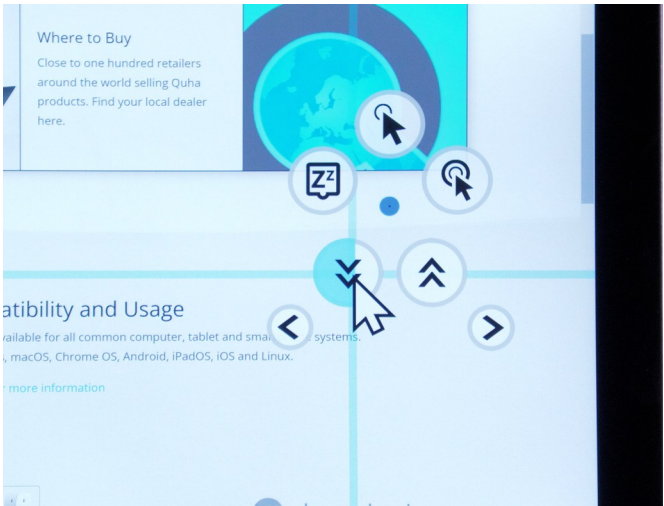

# QUHA PUFO+ - PUFF SWITCH FOR QUHA ZONO AND ZONO 2

Quha Pufo+ is a wireless puff switch for the Quha Zono and Zono 2 gyroscopic mice.

Quha Pufo+ puff switch is quickly connected to Quha Zono — no additional power or configurations are needed. The mouthpiece is easily adjustable to each user and replaceable for hygiene.

P/N	PRODUCT
181	QUHA PUFO+ WITH HEADBAND
182	QUHA PUFO+ WITH HEAD MOUNT KIT



# QUHA KAJO BUTTON - ABILITY SWITCH



Quha Kajo Button is a sensitive and durable Assistive Technology switch designed for everyday use. It can be used in any situation that requires switch access.

Quha Kajo Button has been designed to endure hard use and no compromises have been made with the quality of components. It is mountable using the thread inserts in the bottom.

Quha Kajo Button switches have distinctive auditory and tactile feedback. They are available in three different sizes - 25 mm (1"), 50 mm (2") and 100 mm (4").

# QUHA KAJO SOFT BUTTON - PILLOW SWITCH



Quha Kajo Soft Button is a soft switch with a large soft surface. Because of the soft foam layer under the cover the Soft Button can be safely pressed with head movement.

The safety pin in the bottom of the Soft Button ensures that the button stays where it is supposed to stay during the night. It can be attached to a pillow or other fabrics. The safety pin is detachable.

The cover of the Soft Button is dirt resistant and skin friendly. The cover can be washed (up to 60°C/140°F).

P/N	PRODUCT
250	QUHA KAJO BUTTON 25 MM
251	QUHA KAJO BUTTON 50 MM
253	QUHA KAJO BUTTON 100 MM
252	QUHA KAJO SOFT BUTTON

MORE INFORMATION  
ON ALL QUHA  
PRODUCTS:  
[www.quha.com/products](http://www.quha.com/products)





Designed in Finland by Quha

[www.quha.com](http://www.quha.com)

Copyright 2025 Quha oy.

Quha, the Q-mark, Quha Zono, Quha Zono 2, Quha Zono X, Quha Sento, Quha Pufo+, Quha Dwell and Instant Dwell are trademarks or registered trademarks of Quha oy.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Quha oy is under license. Other trademarks and trade names are those of their respective owners.

iPhone, iPad and AssistiveTouch are trademarks of Apple Inc., registered in the U.S. and other countries. iOS is a trademark or registered trademark of Cisco in the U.S. and other countries and is used under license.